





Fast MAb Separation – Zenix™ 300, 7.8x200mm

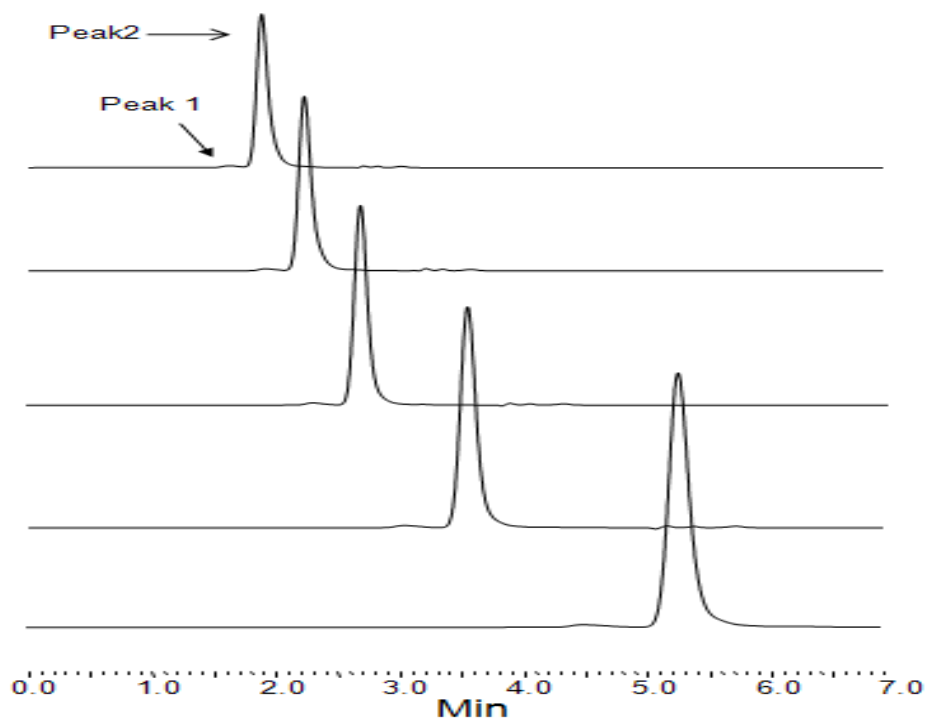

Mobile Phase: 150 mM phosphate buffer, pH 7
Flow Rates: 1.0 mL/min, 1.5 mL/min, 2.0 mL/min, 2.5 mL/min, 3.0 mL/min
Temperature: Ambient
Detection: UV 214 nm
Injection Volume: 10 μ L
Sample: monoclonal antibody

Zoom – in

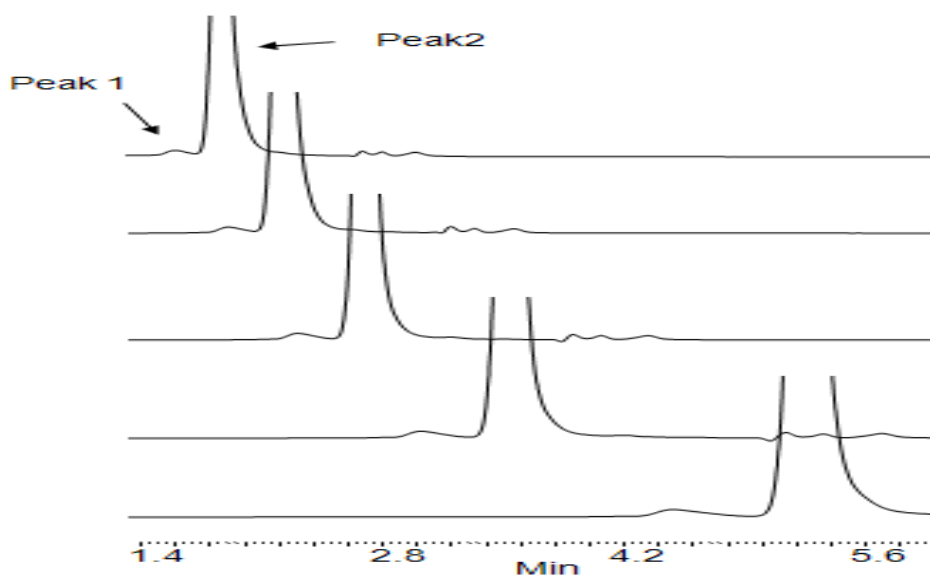

Keywords: Size exclusion, Zenix, monoclonal antibody, fast separation, zoom-in